

VISTA® UPPER SPINE


Raising the Level of Patient Care

Vista Cervical Collar

CODE L0172 APPROVED

The fully adjustable Vista Cervical Collar helps improve patient care while saving time, money and storage space. Built on the design of the trusted Aspen Cervical Collar, the Vista Collar provides motion restriction⁴ while helping to defend against skin breakdown. With its patented height adjustment technology, the right size is always on hand.

4. Comparative Motion Allowed, Vista vs. Aspen, Aspen Medical Products in-house testing, December 2006.

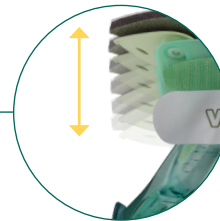


One Size Adjustable

Patented sizing system offers six height adjustment settings to easily accommodate a wide range of anatomies and allows easy adjustments on the patient.



Dial Height Adjustment

Simply pull and turn the Height Adjustment Dial to select one of six height settings.

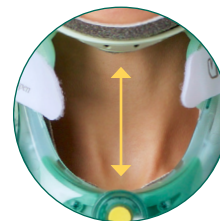

Extra Large Access Window

The extra large tracheal aperture provides functional access for airway management and other essential procedures.



Vista TX Collar

CODE L0174 APPROVED

Encompasses all of the features and benefits of the popular Vista Cervical Collar with an easily detachable thoracic extension, which provides added coverage in the thoracic region.



Vista ICU Back Panel

Pressure Reducing Pad

Unique dual-density foam disperses pressure and helps to reduce skin breakdown: Low density ether foam cradles the patient's head to aid in pressure distribution. Higher density polyurethane foam provides added support and structure.



Extended Posterior Pad

Provides occipital support, isolates from bed shear when the patient is repositioned and can be removed once patient becomes ambulatory.

Vista ICU Back Panel sold separately.

The Aspen Advantage

Aspen Medical Products builds the finest spinal braces on the market. With an unparalleled commitment to Research and Development, Aspen has brought to market numerous unique and innovative features, marked by the Aspen Advantage logo. These features yield real benefits for practitioners and have helped countless patients get back to the activities they care about. Quality products you can count on with features that make a positive difference. That's the Aspen Advantage.



THE ASPEN ADVANTAGE SYMBOL

An innovative feature available only from Aspen

VISTA Upper Spine Braces

Innovative Options for Acute Restriction



Aspen's premier Vista® Upper Spine line utilizes an innovative design and unparalleled technology to enhance patient care. Motion restriction, effective pressure distribution and one size adjustability work together to help improve patient outcomes, minimize inventory needs and aid in eliminating skin breakdown. Based on the degree of motion restriction required, Vista offers a comprehensive line of fully adjustable products to accommodate the patient's clinical conditions and progression of care.

Multi-Award Winning Cervical Care

The unique and innovative design of the Vista Collar and Vista CTO earned both products the prestigious **SPINE TECHNOLOGY AWARD** for Best New Technology in Cervical Care. The Vista Collar also received the TechAmerica **TERMAN AWARD** for the best Medical Device. Aspen is proud to have won these awards and even prouder to be the cervical orthoses of choice in the majority of Trauma Centers across North America, working together to provide superior patient care and improved outcomes.



Aspen Acute Restriction Braces Provide a Controlled Environment for Healing

Aspen Acute Restriction braces provide significant motion restriction by acting as a kinematic restrictor¹, for patients needing support for healing, and can be safely and quickly applied. Gross motion restriction is often desired in treating injured patients recovering from trauma or surgery². The goal is to support healing³, limit recovery time and save money by expediting patient discharge times, accelerating turnaround times and minimizing health risks.

1. Cholewicki J, Lee AS, Reeves NP, Morrisette DC. "Comparison of trunk stiffness provided by different design characteristics of lumbosacral orthoses." *Clin Biomech (Bristol, Avon)*. 2010 Feb; 25(2): 110-114., doi: 10.1016/j.clinbiomech.2009.10.010.
2. Cholewicki, Jacek, et al. "Comparison of Motion Restriction and Trunk Stiffness Provided by Three Thoracolumbosacral Orthoses (TLSOs)." *Journal of Spinal Disorders & Techniques*, vol. 16, no. 5, Oct. 2003, pp. 461-468.
3. Cholewicki, Jacek. "The Effects of Lumbosacral Orthoses on Spine Stability: What Changes in EMG Can Be Expected?" *Journal of Orthopaedic Research*, vol. 22, no. 5, 2004, pp. 1150-1155., doi:10.1016/j.orthres.2004.01.009.

Vista MultiPost Collar

CODE L0180 APPROVED

The Vista MultiPost Collar with pivoting occipital panels adds an extra level of support for greater motion restriction.⁵ Like the original Vista Collar, the Vista MultiPost Collar is one size adjustable, virtually eliminating waste associated with collar sizing errors. It also features Aspen's proven cotton-lined pads designed to enhance skin care by reducing patient contact points in the occipital area.



Pivoting Occipital Panels Reduce Pressure [Ⓐ]

Padded panels pivot bilaterally and self-adjust to cradle almost all head shapes for increased support and to help maximize patient comfort. Reduced patient contact on occipital midline helps limit pressure points and heat buildup in areas prone to skin break down.




Increased Visibility for Improved Care

Innovative design permits greater visibility, enhancing the ability to inspect the skin without removing the collar.


Easily Steps Up to Support Patients' Condition [Ⓐ]

If additional motion restriction is required, both the Vista CTO and Vista CTO4 Upgrade Kits attach easily to the Vista MultiPost Collar.

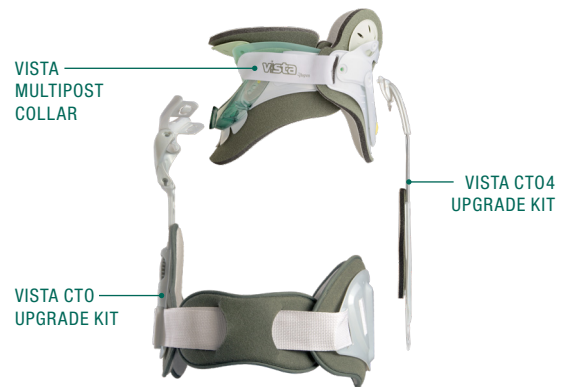

Promotes Effective Skin Care [Ⓐ]

Cotton-lined Aspen pads wick moisture from patient's skin while clickable polyurethane foam helps ensure optimal pressure distribution.⁶

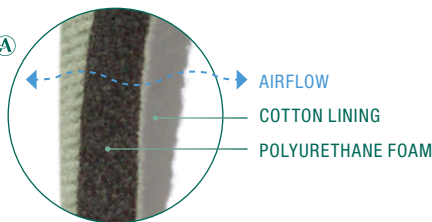

Fully Adjustable for the Perfect Fit

Height adjustment technology in both the front and back allows clinicians to effectively fit individual anatomies and multiple conditions.



Vista Collar Front Panel [Ⓐ]

Features the Vista Collar front panel with the patented, easy height adjustment system for one size adjustability and an extra large tracheal aperture for improved access and airway management.

5. Ivancic, PC., et al. Do Cervical Collars and Cervicothoracic Orthoses Effectively Stabilize the Injured Cervical Spine? Spine, 2013 Jun 1;38(13):E767-74. *Vista TS & Vista TS4 now Vista CTO & Vista CTO4.
6. Blaylock, B., et al. Solving the problem of pressure ulcers resulting from cervical collars: a step-by-step approach. Ostomy Wound Management. 1996 May; 42(4): 26-33.

Vista CTO and Vista CTO4

Designed to provide motion restriction throughout the cervico-thoracic spine, the Vista CTO and Vista CTO4 provide significant motion control in all three planes of motion—flexion/extension, lateral bending and axial rotation.⁵ These one size adjustable braces feature multiple adjustment points to help ensure a comfortable fit and the ability to step down to support patient recovery.



Locking Sternal Adjustment [Ⓐ]

Angle of support strut is easily adjusted to fit patient's anatomy then locks into place.

Padded for Comfort

Fully padded at all patient contact points to help maximum comfort, encouraging improved patient compliance.



Flexible Thoracic Adjustment [Ⓐ]

One size fits essentially all patient anatomies and unlocks to adjust during original fit and repositioning.



The Vista CTO4 [Ⓐ]

For enhanced motion restriction, the Vista CTO4 provides increased motion control in all three planes of motion.⁵ The adjustable, low-profile back strut is designed for maximum patient comfort, even in the supine position.



VISTA CTO4



VISTA CTO



VISTA
MULTIPOST
COLLAR

Easily Steps Down to Support Progression of Care [Ⓐ]

As a patient's condition improves, the Vista CTO4 can be stepped down to the Vista CTO then down again to the Vista MultiPost Collar, all without the use of tools.

Product Information

Vista MultiPost Collar	Vista MultiPost Collar	984200
	Vista MultiPost Collar Set with Replacement Pads	984202
Vista TX Collar	Vista TX Collar	984100
	Vista TX Collar Set with Replacement Pads	984102
Vista Cervical Collar	Vista Collar	984000
	Vista Collar Set with Replacement Pads	984002
	Standard Back Panel	984006
	Large Back Panel	984008
Vista ICU Back Panel	Vista ICU Back Panel	984600
Vista CTO4 · 4 Post CTO	Vista CTO4	984550
Vista CTO · 2 Post CTO	Vista CTO	984510
Upgrade Kit	Vista CTO4 Upgrade Kit · Posterior	984554
	Vista CTO Upgrade Kit · Anterior	984504
Replacement Pads	Vista MultiPost Collar Replacement Pads	984220
	Vista TX Collar Replacement Pads	984120
	Vista Collar Replacement Pads	984020
	Vista Collar Large Back Panel Replacement Pad	984028
	Vista ICU Back Panel Replacement Pad	984620
	Vista CTO4 / Vista CTO Replacement Pads	984520